

WHEELCHAIR TORSO SUPPORT SELF-RELEASE POSITIONING HARNESS

Model #: WCH-1 • Made in PRC

SINGLE PATIENT USE

IMPORTANT INFORMATION

Please read all instructions, warnings and precautions before use. Correct application is essential for proper functioning of the product.

Use only on the person it was provided to by your facility and only for the intended use.

INTENDED USE

The Secure[®] Torso Support Self-Release Wheelchair Positioning Harness is intended to be used to help keep your resident or patient properly positioned in wheelchair, prevent forward sliding and possible attendant falls and getting up unassisted.

WARNINGS AND PRECAUTIONS


- Do not use on a resident or patient who is or may become highly aggressive or agitated.
- NEVER use on a resident or patient if an I.V. or wound site could be compromised by the harness.
- NEVER use as a restraint.
- Always secure waist strap to a portion of the wheelchair that moves with the patient.
- Must be worn over a resident or patient's clothing.
- Check the resident or patient regularly to ensure circulation is not compromised.
- ALWAYS ensure the resident or patient can easily self-release the two shoulder straps and waist strap. If the user cannot easily remove the complete harness DO NOT USE.

INSTRUCTIONS FOR USE

1. Wrap the waist strap around the resident or patient and the wheelchair. NOTE: Resident or patient should be sitting up while device is being applied. The waist belt closure MUST be position in the front of the resident or patient so the waist belt can be easily opened and released by the user.
2. Secure hook and loop closure.
3. Cross the shoulder straps in an X behind the wheelchair back.

4. Bring the straps over the resident or patient's shoulders and securely attach straps to the front of the belt. NOTE: Secure excess slack in straps.

5. Do not over tighten. Two fingers should fit between the device and the patient to ensure circulation is not compromised.

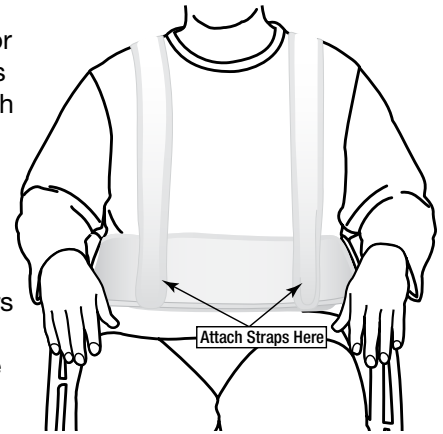

SIZING

Universal

CLEANING AND/OR MAINTENANCE

Close any hook and loop fasteners before washing to help prevent lint build up and pilling. For maximum product life and to help prevent any straps from tangling, wash in a mesh laundry bag. Hand or machine-wash with a mild fabric detergent in warm water up to 60°C / 140°F. Bleach as directed on container.

Do not iron. Lay flat and open and allow to air dry, or tumble dry on low heat for up to 40 minutes. Inspect the harness after drying to ensure it is intact. Hook and loop fasteners may collect lint after repeated use or laundering, reducing grip strength. As needed, use a stiff bristle brush to remove lint from the "hook" side. Test hook and loop adhesion before each use and discard product if hook-and-loop does not hold securely. If the product is damaged, do not use it and replace it with a new one.

SIX MONTH LIMITED WARRANTY

Personal Safety Corporation warrants this product to be free from factory defects and workmanship for a period of six months from date of purchase. Warranty is void if this product is misused, abused or damaged. No warranty or responsibility of use expressed or implied.



Please see our complete line of Secure[®] fall/wandering prevention & patient safety products at www.SecureSafetySolutions.com

Highest Quality...Lowest Prices...Guaranteed!

Personal Safety Corporation • Hiawatha, Iowa 52233 • 800.373.2873