

# 14 Day P.A.D.S. Chair Sensor


## The Personal Accident Deterrent System

Can be used with Secure<sup>®</sup> Universal Alarm (SUA-1 or SUA-1Y),  
Secure<sup>®</sup> Code Alarm (SCA-1), Secure<sup>®</sup> 120 Alarm (SUA-120),  
Secure<sup>®</sup> Fall Alarm Monitor (CAS-FAM)

Model No. CHAIRPAD-14 (12" x 12") • Made in PRC

**CAUTION!**  
ALWAYS remove sensor pad pin from alarm by pulling straight up on the white "PULL" tab. DO NOT attempt to remove pin any other way or you may damage pin and void the warranty. Thank you.

ALARM AND HOLDER MUST BE AFFIXED OUT OF REACH OF RESIDENT  
  
U.S. PATENT NO. 6,847,301

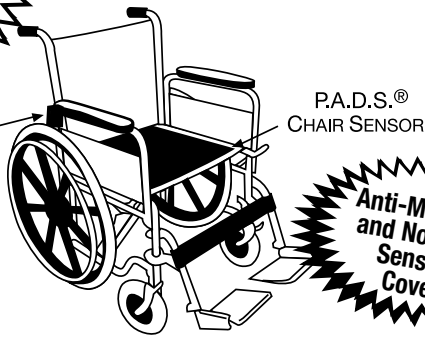

### Instructions

1. The alarm and sensor pad should be checked prior to each use for proper functioning. The "Test/Reset" button also functions as a battery tester.
2. Place P.A.D.S.<sup>®</sup> Chair Sensor in center of chair seat directly under buttocks of user. If chair cushion is used sensor pad must be placed on top of chair cushion directly under buttocks of user.
3. Install batteries (included) in Alarm by opening cover with screwdriver. Put sensor pad connecting pin into opening at top, front or pad jack opening of alarm. **CAUTION! ALWAYS remove sensor pad pin from alarm by pulling straight up on the white "PULL" tab. DO NOT attempt to remove pin any other way or you may damage pin and void warranty.**
4. Affix Alarm Holder securely to bottom or back of chair out of sight and reach of the user. Simply wrap Velcro<sup>®</sup> strap with rubberized non-slip backing around railing or frame of chair. Holder placement should be easily accessible to staff. Put alarm in holder and fasten Velcro<sup>®</sup> closure of holder to Secure<sup>®</sup> alarm in holder. Secure excess cord to avoid damaging.
5. Carefully assist user into chair.

6. When using the Secure<sup>®</sup> Universal Alarm (SUA-1 or SUA-1Y) with this pad push the "Test/Reset" button on alarm to test. The alarm will sound 4 times and shut off to indicate proper functioning. Alarm will allow for normal movement. However, if user gets off sensor pad the alarm will sound to alert staff. When user gets back on pad the alarm will automatically reset. **PLEASE NOTE: If using a Secure<sup>®</sup> brand alarm other than the Secure<sup>®</sup> Universal Alarm (SUA-1 or SUA-1Y) with this pad please follow instructions for the particular alarm. Instructions for all Secure<sup>®</sup> brand alarms can be printed from our web site at [www.SecureSafetySolutions.com](http://www.SecureSafetySolutions.com).**
7. When staff removes user from pad for whatever reason the alarm will sound. To silence alarm press "Test/Reset" button per alarm's instructions.

### Cleaning Instructions

Use the same cleaning method as you use to clean vinyl mattress covers.

### Water and Incontinent Proof

The P.A.D.S.<sup>®</sup> Sensor has been specially designed to be water and incontinent proof. If sensor is exposed to extreme heat, punctured or torn it will no longer be water and incontinent proof and should no longer be used. Warranty will be voided. Never submerge sensor and pad completely in water.

### 14 Day Limited Warranty

**Personal Safety Corporation warrants this product to be free from factory defects in materials and workmanship for a period of 14 days from date below. Void if this product is misused, abused or punctured. No warranty or responsibility of use expressed or implied.**

## CAUTION!

This product is designed as a monitor, alerting the staff when a patient has risen from their bed or chair. This product is not a restraint and therefore will not prevent a patient from getting out of their bed or chair. This product is to be used in conjunction with a total fall prevention program.

## DO NOT ASSUME!

This product must be tested each time it is used with a patient to make sure it functions properly.



Please see our complete line of Secure<sup>®</sup> fall/wandering prevention & patient safety products at [www.SecureSafetySolutions.com](http://www.SecureSafetySolutions.com)

**Highest Quality...Lowest Prices...Guaranteed!**

**Personal Safety Corporation • Hiawatha, Iowa 52233 • 800.373.2873**